

C1QBP Polyclonal Antibody

Catalog No: #40661

Package Size: #40661-1 50ul #40661-2 100ul

Orders: order@signalwayantibody.comSupport: tech@signalwayantibody.com

Description


Product Name	C1QBP Polyclonal Antibody
Host Species	Rabbit
Clonality	Polyclonal
Purification	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.
Applications	WB IHC
Species Reactivity	Hu Ms Rt
Specificity	C1QBP Polyclonal Antibody detects endogenous levels of C1QBP protein.
Immunogen Description	Synthesized peptide derived from the Internal region of human C1QBP.
Target Name	C1QBP
Other Names	C1QBP; GC1QBP; HABP1; SF2P32; Complement component 1 Q subcomponent-binding protein; mitochondrial; GC1q-R protein; Glycoprotein gC1qBP; C1qBP; Hyaluronan-binding protein 1; Mitochondrial matrix protein p32; p33
Accession No.	Swiss-Prot: Q07021NCBI Gene ID: 708
Uniprot	Q07021
GeneID	708;
SDS-PAGE MW	32kd
Concentration	1mg/ml
Formulation	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide.
Storage	Store at -20°C/1 year

Application Details


Western Blot: 1/500 - 1/2000.

Immunohistochemistry: 1/100 - 1/300.

Images

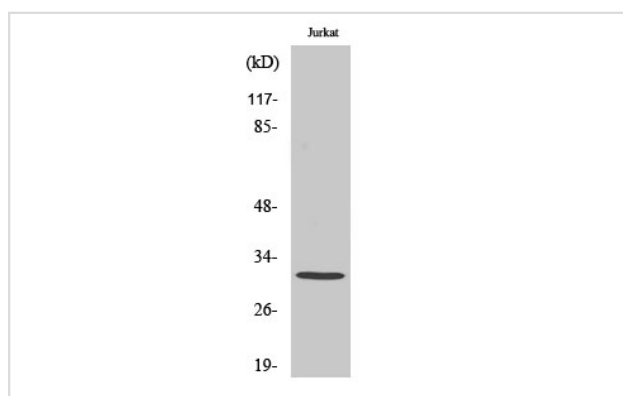

Western Blot analysis of MCF7 cells using C1QBP Polyclonal Antibody

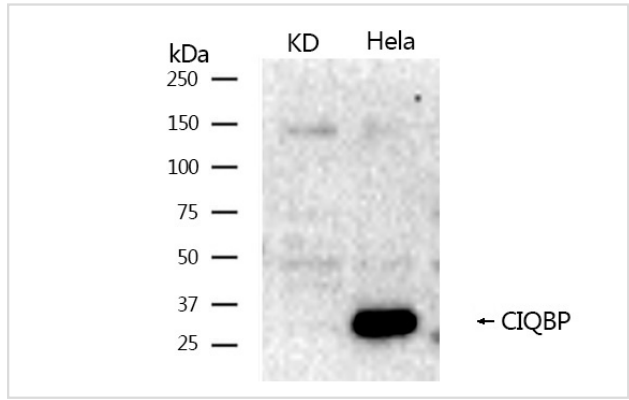

Western blotting analysis using C1QBP Polyclonal Antibody #40661.

Note: This product is for in vitro research use only